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## EN - Clipp® user instructions

**PLEASE NOTE** the Clipp® weaning device must not be used on any child's dummy other than the one provided in the kit and certified as compatible.
**PLEASE READ THE INSTRUCTIONS CAREFULLY**


CLIPP® is a CE medical device for Orthodontic prevention from 2 years of age, without BPA (bisphenol A), consisting of a dummy and 5 identical shields; all of these 6 elements are intended to completely wean the child off the dummy in 30 days (minimum). It is based on the principle of gentle weaning in stages, at the child's own pace and with no going back. **Please follow the instructions below closely.**

**The CLIPP® device consists of:** 1 dummy, 5 weaning shields (identical) and detailed instructions.

**Therapeutic indications:** this device is intended for children who show signs of dependency on the dummy or whose parents wish to stop it from being used. It can be used from 2 years of age, is recommended in the ideal period of around 3 years of age and is strongly recommended in cases of late non-nutritive sucking after 3 and a half years of age.

**Contraindication:** The Clipp dummy should not be chewed. If the child chews his/her dummy, there is a significant risk of damage to the device, which is not suited for this use (intended for late sucking). We advise offering the child objects that are more suitable for chewing or an orno-my functional assessment in order to find the origin of the problem.

**Instructions for use:** For each of the following stages you can double the recommended time if necessary depending on your child's pace or just slow down at certain stages whilst adhering to an overall maximum of 60 days.

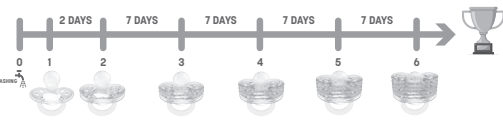**PLEASE NOTE ALL STAGES (1 to 6) MUST BE MANAGED AND CONTROLLED BY AN ADULT.**

**For the safety of the child:** once they are clipped on, the shields can no longer be removed.

1-Wash the 6 elements in hot water before use: the dummy and the 5 weaning shields, (all identical)

Replace the child's dummy(ies) with the Clipp® dummy(ies) only WITHOUT the weaning device (without shield) then wait 2 to 3 days before moving on to stage 2 below.

At this point the dummy is similar to other dummies on the market but it is imperative to remove all other dummies available to the child. Stopping all sources of sucking, such as the bottle, comforter or other items is also recommended, in order not to hinder the weaning process.

2-Add the first shield (no particular orientation for this first one but there will be a specific methodology starting from the next one)

It is important to ensure that THE 3 CLIPS of the first shield are securely locked onto the base dummy. The dummy + shield 1 form a non-separable unit whose sucking teat is approximately 4 mm shorter. This unit will enable the weaning phase to begin. Wait 7 days then proceed to stage 3 below.

3-Add the second shield PLEASE NOTE: THERE IS A SPECIFIC WAY OF DOING THIS! See the graphic on the reverse side or the linked video: QR code.

It is important to securely lock the 3 clips of the second shield into the 3 ROUND available holes of the previous shield. The dummy + shields 1 and 2 form a non-separable unit whose sucking teat is a little shorter. Wait for a minimum of 7 days then proceed to stage 4 below.

4-Add the next shield using the same principle, at the child's own pace. PLEASE NOTE THE METHODOLOGY see the graphic on the reverse side / It is important to securely lock the 3 clips of this new shield into the 3 ROUND available holes in the previous shield already in place. The dummy + the shields form a new non-separable unit reducing the sucking teat still further. Wait for a minimum of 7 days then proceed to stage 5 below (take your time, the next stage is more difficult for the child, it is important to proceed at the child's own pace).

5-Repeat stage 4 for the next shield. At this stage, it becomes very complicated for the child to keep the Clipp® dummy in the mouth, and the child gradually abandons it on his/her own. Wait 7 days then proceed to stage 6 below.

6-Fix the last shield (fifth) using the same principle. The child can no longer use the dummy, it is "locked", but he may keep it for sleeping in the same way as a soft toy or decide to throw it away to officially mark the definitive end.

Congratulations! The weaning has finished, it's a success!

**If using several Clipp® devices simultaneously: stages 1 to 6 of the instructions should be done at the same time on all the dummies.**

From an orthodontic point of view as well as for this weaning dummy we advise trying to limit the use of the dummy to specific times such as at bedtime, during a nap or during periods when the child is very agitated, for example. It should not be used systematically; completion of the weaning process depends on this.

In the event of a gag reflex before 3 years of age, the first shield can be added immediately to reduce the length of the sucking teat (stages 1 + 2). If your child has a tendency to chew and damage dummies, he risks doing the same with the Clipp® dummy; you can try and train the child and provide a second kit to be able to complete the weaning.

**PLEASE NOTE!** Before each use, carefully examine the dummy and pull it in all directions. Discard it if there is the slightest sign of deterioration or fragility. An adult should regularly check that the sucking teat is not damaged and that the clips which hold the shields in place are not broken. The shields should not be removable once correctly installed.

If one day you find that the child has the whole dummy in his/her mouth, DO NOT PANIC - the dummy cannot be swallowed, it is manufactured so as not to be a danger to your child. Remove the dummy very carefully.

Do not wash in the dishwasher, do not disinfect in a pressure cooker - there is a risk of material damage from detergents or from uneven heating. Do not expose a dummy directly to sunlight or to heat and do not leave it in a disinfectant product ("sterilisation solution") for longer than the recommended time as this could weaken the dummy. Do not soak in sugary food products or medical solutions - there is a risk of cavities and/ or material damage!

**Wash the device regularly with hot water throughout the weaning phase; you can use a small bottle brush to clean between the shields.**

-----

**For your child's safety - WARNINGS!**

The product must not be attached with anything around the neck; never attach it to cord, ribbon, string or items of loose clothing.

Only if necessary, use a dummy clip; ensure that it complies with the regulations currently in force (standard EN 12586), 22 cm maximum, never extend the dummy clip!

Avoid attaching a soft toy to the dummy in order to limit the forces exerted on the teeth and so limit malocclusions.

-----

CE product - all 6 elements are certified 0% BPA (bisphenol A) in compliance with current regulations - from 2 years of age.

The dummy alone (main element) is certified compliant: BS EN 1400/2013+A2:2018 Child use and care articles - Soothers for babies and young children - Safety requirements and test methods. BS EN 71-3:2019+A1:2021 Safety of toys - Part 3: Migration of certain elements.

REGULATION (EC) No. 1907/2006 (REACH) of the European Parliament and of the Council of 18 December 2006, including its amendments.

device distributed by: Biomimétik 2 bis rue des Marettes 28320 Ymeray instruction notice v7 11.2023

Recommendation of the German Federal Institute for Risk Assessment - Recommendation XV, silicone

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